




# Air Curling

332-C1904

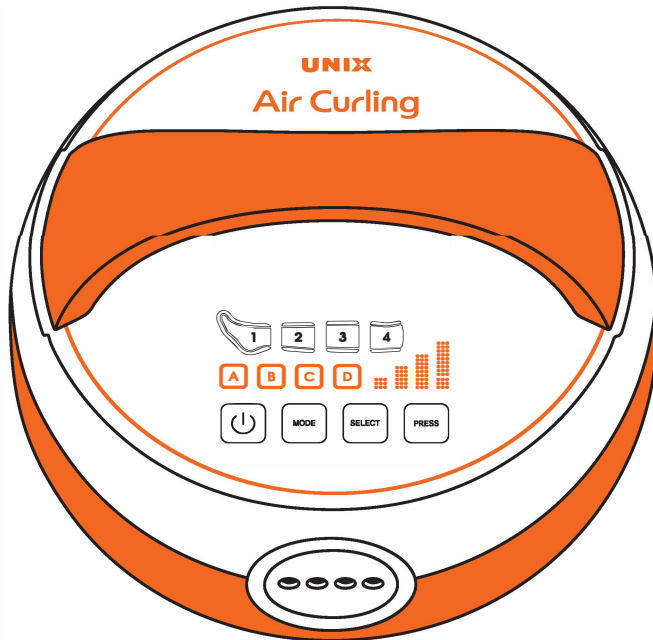

---

Read this manual carefully before use and save it for later use.  
This device is NOT FOR MEDICAL USE.

# CONTENTS



## Air Curling

332-C1904

<b>CONTENTS</b>	02
<b>Cautions on safety</b>	03
<b>Warnings on use</b>	04
<b>Cautions on use</b>	05
<b>Prohibitions on use</b>	06
<b>Product &amp; Accessories</b>	07
<b>Component Names</b>	08
<b>Display Window</b>	09
<b>Function buttons</b>	10~11
<b>How to use this product</b>	12~14
<b>Maintenance</b>	15
<b>Specifications &amp; Warranty Card</b>	16
<b>Troubleshootings</b>	17
<b>MEMO</b>	18~19



Air compressed cuff with four chambers



LED Display



Various massage programs



11ℓ/1min Air pump



4 Levels of pressure adjustment



Wireless Control

# Cautions on safety

4 signs in this instruction manual are adopted to prevent any injuries or damages and operate it properly using this product.  
Please read thoroughly and keep it in a near place to read it next time.

-  **WARNINGS**  
.....
-  **CAUTIONS**  
.....
-  **PROHIBITIONS**  
.....
-  **OBSERVANCE**  
.....

## USE PURPOSE

It relaxes your legs, arms and waist by inflating and deflating the air cuffs.

# Warnings on use



## WARNINGS



- 01) Consult with your doctor if you are under a medical treatment.
- 02) Do not take it apart.(It can cause a problem or electrical shock.)
- 03) It is recommended to use 3 times a day and one time is for 15 minutes. Do not exceed the recommended time.
- 04) It can hurt yourself if you use it continuously for a long time or use it at a strong power level at initial operations.
- 05) Consult with your doctor if you felt uncomfortable with yourself.
- 06) In case of power interruption, Turn off and plug out.
- 07) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- 08) If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- 09) Use the appliance only for its intended use as described in the manual.
- 10) Do not plug in or pull out with wet hands.(It may cause an electric shock.)
- 11) Do not use this product in a wet place or hot places. (It may cause a defect, electric shock, fire and etc.)
- 12) Only operates this product after you wear the air cuffs. (The air cuff can be damaged.)

# Cautions on use



CAUTIONS



OBSERVANCE



- 01) Neither give a strong impact on the product nor drop it. It may cause a damage.
- 02) Keep the product close to you to operate it easily.
- 03) Open the air boots and be careful that nothing is inserted in the zipper.
- 04) Do not move or walk while wearing the air boots.
- 05) If any trouble or defect is found, Stop using the product immediately and turn it off and contact your local distributor.
- 06) Be careful that water or dirty materials do not go into the cuffs.
- 07) After use, The air in the air boots goes out fully and then keep it in a proper place.
- 08) Do not put other things on this product.
- 09) Plug it off in case of not using for a long time.
- 10) Be careful that the air hoses are not twisted or pressed during use or keeping.
- 11) Keep this product in a proper place.
- 12) Keep the cuffs away from needles or sharp objects. (It may cause a fatal damage.)
- 13) Please cover the zipper area with a cuff cover inside and pull the tab.
- 14) Call your local distributor for repair.

# Prohibitions on use



## WARNINGS



## PROHIBITIONS

“Please consult with a doctor before use if you have and had the following diseases.”

It is prohibited to use this product if you use a medical device.

**People with the following disease or condition are prohibited to used this product.**

- |                                       |   |
|---------------------------------------|---|
| 1) Acute disease, High blood pressure | 9) Cerebral hemorrhage  |
| 2) Diabetes                           | 10) Tuberculosis  |
| 3) Hemophilia                         | 11) Drinking alcohol or taking a medicine   |
| 4) Atopic dermatitis                  | 12) The pregnant  |
| 5) Shedding of blood                  | 13) Malignant tumor   |
| 6) Menstruation                       | 14) Under bad health condition  |
| 7) Heart disease                      | 15) Heavy fatigue   |
| 8) Allergy                            | 16) People who have been diagnosed by doctors<br>as inappropriate to use this machine |

### **Other prohibited things**


- 1) Take things out from your clothes before use.
- 2) Consult with a doctor if you feel something wrong with your body.
- 3) Consult with a doctor if you feel itchy or have a rash on your body.
- 4) Do not use this product after putting medicine.
- 5) Keep out of children's touch.
- 6) Do not use the product for other purposes.

# Product & Accessories


## Device & Accessories

Device / 1 pair of leg cuff / Remote controller / Foot plate


Device

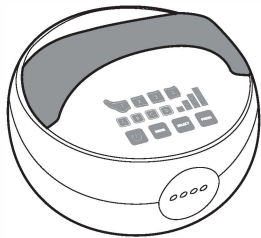

1 pair of leg cuff

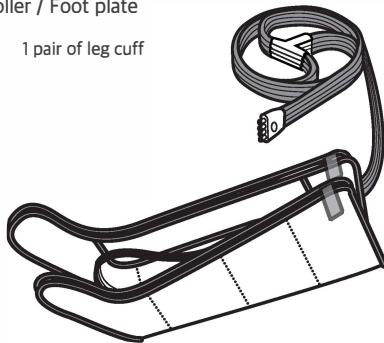

Foot plate

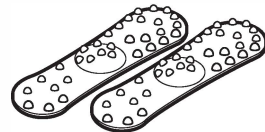


Remote controller



---

## Option


Arm Cuff / Waist cuff / Half Pants cuff / Extender

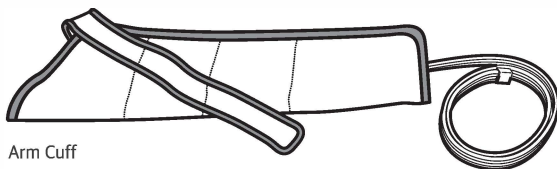


Arm Cuff

Waist Cuff


Extender

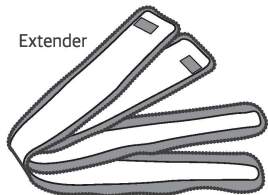

Half Pants Cuff

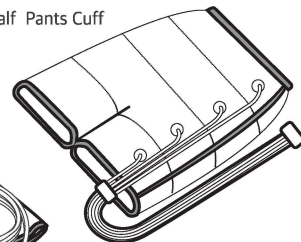

# Component Names

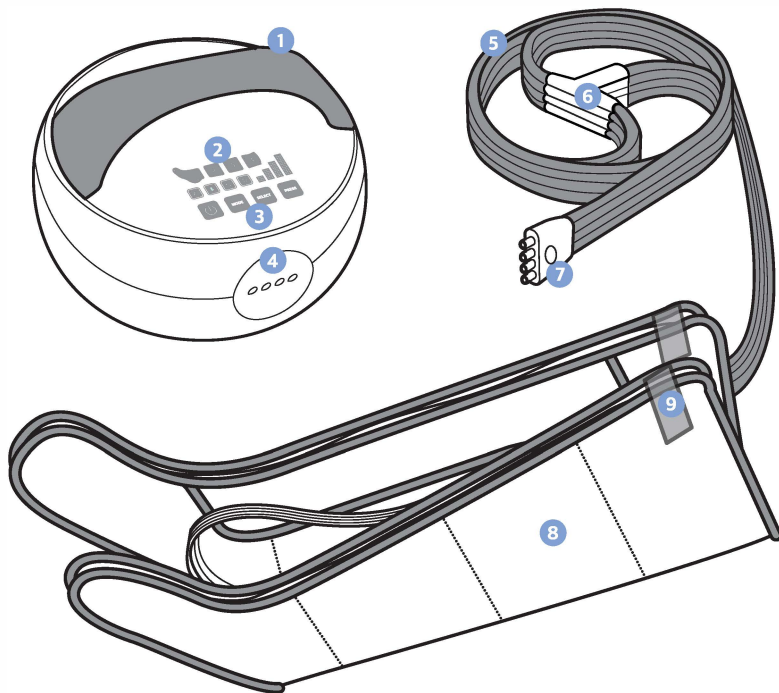


- 1 Portable grip
- 2 Operation display window
- 3 Operation buttons

- 4 Air socket
- 5 Silicone hose
- 6 T connector

- 7 Hose plug
- 8 Leg cuff
- 9 Velcro wrap



# Display & Control Panel

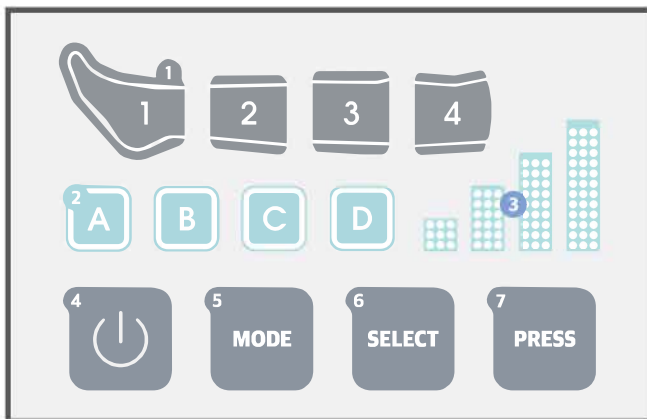


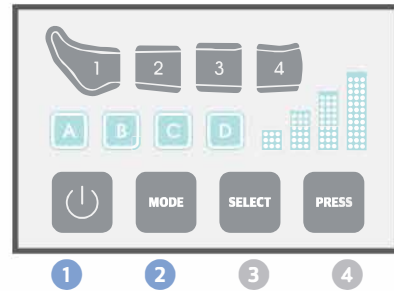
Push the switch gently

- 1 Compression area display**  
It displays the status of inflation and deflation by LED light.  
ON: Inflation or Keeping the inflation  
OFF: Inflation or Deflation
- 2 Display of MODE A, B, C, D**  
MODE A is selected when you turn it on, and you can choose another mode.
- 3 Pressure level Display**  
It displays the status of pressure level.
- 4 ON/OFF Button**  
It displays the status of ON and OFF.
- 5 MODE Button**  
4 Modes can be chosen.
- 6 SELECT Button**  
It can select desired chamber you want.
- 7 Pressure Adjustment Button**  
It can adjust pressure by 4 steps.

# Function Buttons 1


## 1 ON/OFF

- It displays the status of ON and OFF.
- MODE A is selected when you press  button at first and the pressure is Level 1.
- The device will be automatically off after 20 minutes.




## 2 **MODE**

- Press the “MODE” button if you want to change it to B, C and D. (Refer to the program below)
- If you press the “SELECT” button during use, The “MODE” lamp gets off and the lamp on the chambers gets on.

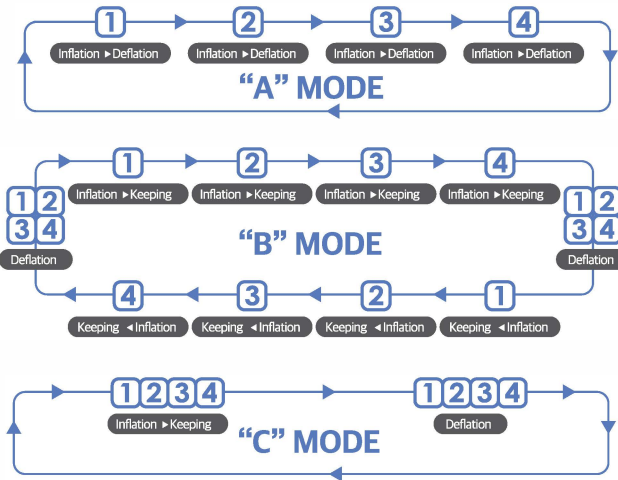


**“D” MODE :** It repeats the modes of A, B and C in order.

### Arm Cuff(Optional)

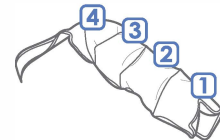


### Waist Cuff(Optional)

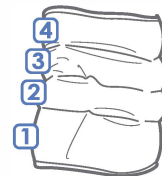

Optional cuffs also have 4 chambers as the same as the leg cuff.


# Function Buttons 2

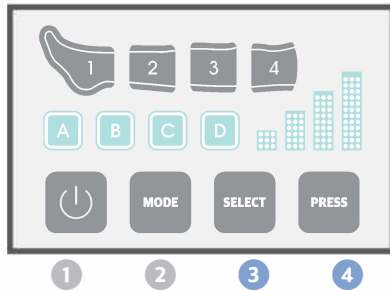

## 3 SELECT

- A user select the chambers to be inflated.
- It repeats inflation and deflation.

[SELECT] The pressure area is changed whenever the “SELECT” button is pressed.



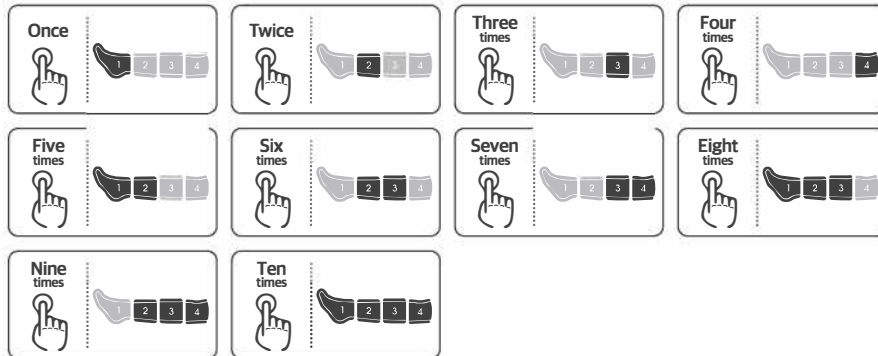


## 4 PRESSURE

- It consists of 4 pressure levels and you can select it as you press the button.
- Deflating function is working temporarily in case the PRESSURE button is pressed.

# How to use this product 1



## 1 Connect the Hose Plug to the Control Unit

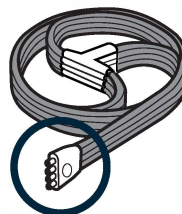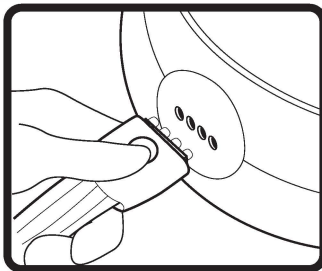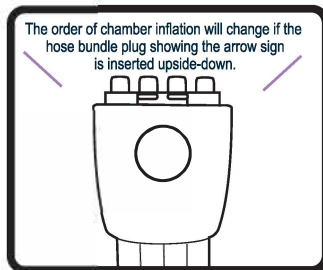

Hose plug

- 1) The hoses should not be twisted or folded.
- 2) Connect the hose plug to the hose plug receptacle on the control unit making sure the arrow sign is facing up.
- 3) The hose plug should be connected fully to prevent air from leaking out.

## 2 Wear a cuff

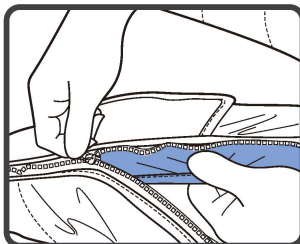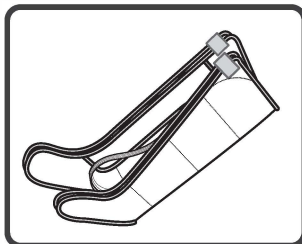


Protecting cover

Cover the protecting cover to protect your skin.

- 1) Wear a light clothes to use this product.
- 2) Check left and Right cuff for leg cuff.
- 3) Pull the zipper up to the end and be careful that your clothes, stocking and others are not caught in zipper.

# How to use this product 2

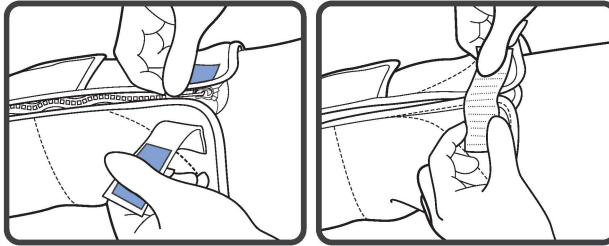

## Velcro

Attach the velcro to prevent any harms on cuffs and your body.

4) Place the zipper flat as in the picture and attach the velcro firmly


### 3 Operate the Product/Remote Control

- 1) Read Page 9 to page 11 in detail to know about the operation and function buttons.
- 2) Place the product in a flat area
- 3) The changed operation starts after deflating the air in the air cuffs for 3 second in case other function buttons are pressed during use




# How to use this product 3

## 1 Connecting an extender(optional)

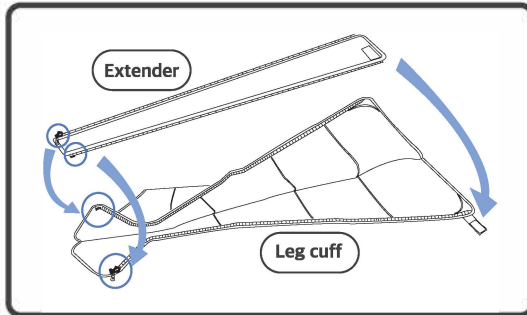


1) Pull the zipper up to the end and attach the velcro firmly.

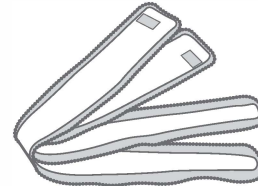

**Extender**

It is recommended for a big person.


## 2 Wearing a waist cuff(optional) / an arm cuff(optional) / Half Pants Cuff(Optional)

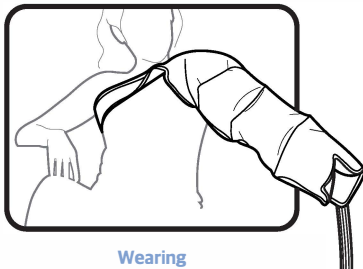

**Wearing  
an arm cuff**



**Wearing  
a waist cuff**

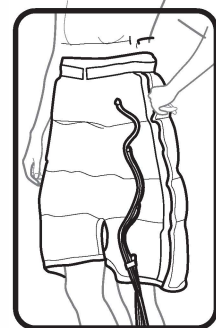

**Wearing  
a Half Pants Cuff**

### **WARNINGS**

- 1) Select the compression strength from low level to high level.
- 2) Do not move around wearing a cuff.

# Maintenance

## 1 After use

- Press ON/OFF button to stop using the product. (Except for automatic turn off after 20 minutes.)
- Disconnect the air hose plug connected to the air socket of the product body.
- Take off the air boots after taking the remained air in the air boots out.
- Disconnect the power cord from an electric outlet.

## The maintenance of a product body

- ### 2
- Keep a product body clean in a cool place.
  - Keep it in a safe place without the danger of drop or shock.
  - Do not clean it with flammable chemicals.
  - Just clean a product body with a soft cloth.
  - Keep it dry after cleaning.
  - It should be stored out of children's touch.
  - Disconnect the power cord if you don't use it for a long time.

## The Maintenance of a cuff

- ### 3
- Clean sweat or water in a cuff with a dry cloth and dry it fully. (Clean it with a soft cloth and do not clean with flammable chemicals)
  - Keep them in a cool place.
  - Do not operate the product without wearing a cuff.
  - Do not keep it close to sharp objects and hot or cold places.
  - Be sure that the connected air hoses are not twisted or folded.
  - Do not try to fold them hard.

# Specifications / Warranty Card

## SPECIFICATIONS

<b>Product Name</b>	AIR LEG MASSAGER
<b>Model Name</b>	332-C1904
<b>Voltage</b>	AC 100~240 V, 50/60 Hz
<b>Power Consumption</b>	30 W
<b>Weight &amp; Packing</b>	2.3 kg, 1set
<b>Product body size</b>	230(W) x 230(D) x 150(H)
<b>Origin</b>	mm MADE IN KOREA

## WARRANTY CARD

<b>Product Name</b>	AIR LEG MASSAGER
<b>Model Name</b>	332-C1904
<b>Date of Purchase</b>	
<b>Warranty</b>	<b>1 year from date of purchase</b>

## CUSTOMER

<b>Name</b>	
<b>Tel.</b>	
<b>Dealer</b>	





# Troubleshooting

## Check Points

Please check the below things before asking for repair service.

Then plug out the power cord in case it is broken.

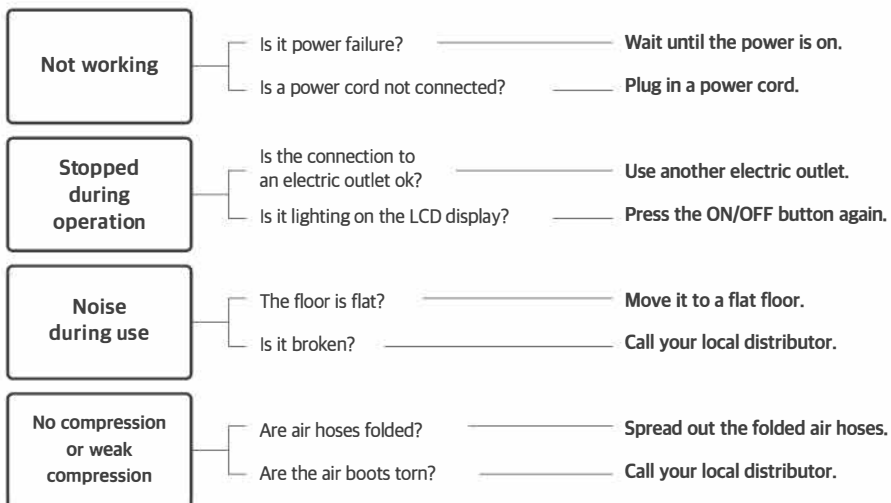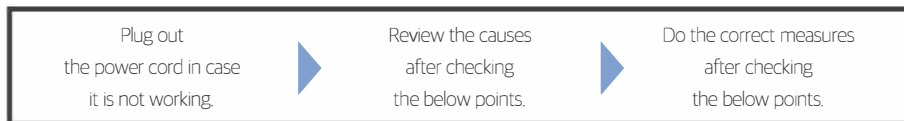

MEMO



A series of horizontal dotted lines spanning the width of the page, providing a guide for writing.

MEMO



A series of horizontal dotted lines for writing, spanning the width of the page.



# Air Curling

332-C1904

Chiroform A/S, Mariendalsvej 22 - 8800 Viborg, Denmark.  
Tel. +45 8661 3611 \* [chiroform@chiroform.dk](mailto:chiroform@chiroform.dk) \* [www.chiroform.dk](http://www.chiroform.dk)

