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
Manual

- DK** Mini Cooper Pro MassageGUN
- EN** Cooper Vibration Massage Device

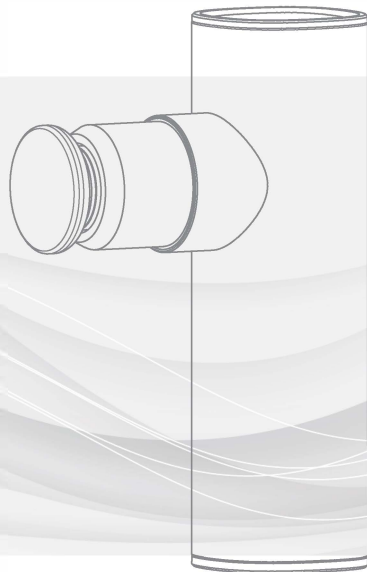● Specifikationer

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
● Pakkeliste

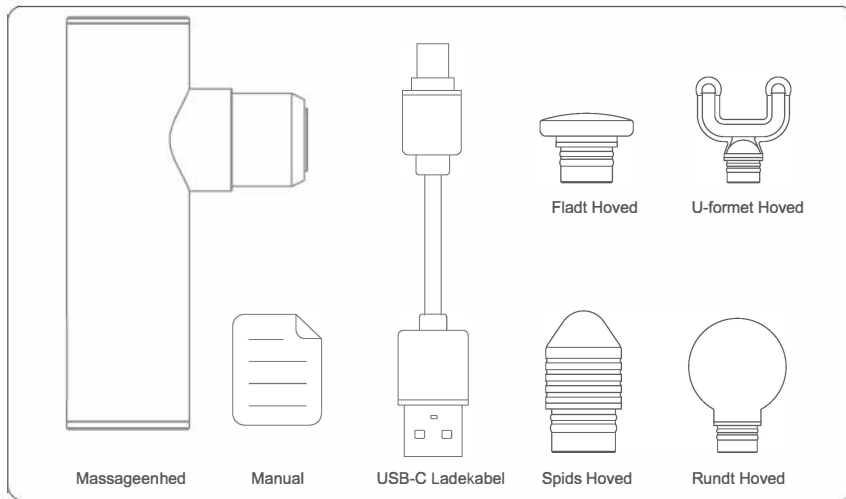Product parameter

Product Name: Vibration Massage Device

Product Model: Cooper

Rated Power: 16W

Peak Power: 60W

Rated Voltage: 8V

Rated Current: 2A

Rated Frequency: 50/60Hz

Output Voltage: 5V

Output Current: 2A

Product Weight: 390g/0.86lbs


The Highest Speed: 3200r/min

Battery Capacity: 2000mAh


Product Size: 75x38x130mm/2.95"x 1.5"x5.12"

(Massage head attachments excluded)

Power Supply: Chargeable Lithium Battery

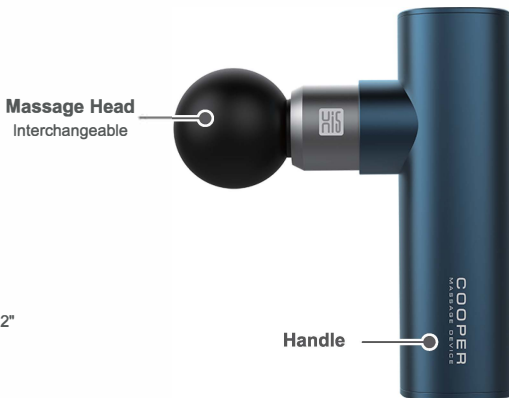Massage Head

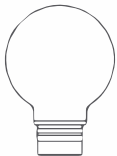Round Head

For full body massage, such as the palm and sole.


Flat Head

For full body relaxation and shaping.

Bullet Head


For deep percussion to meridians and joints.

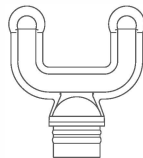U-shape Head

For neck, shoulder and tendon.

🕒 The Way to Change the Massage Head

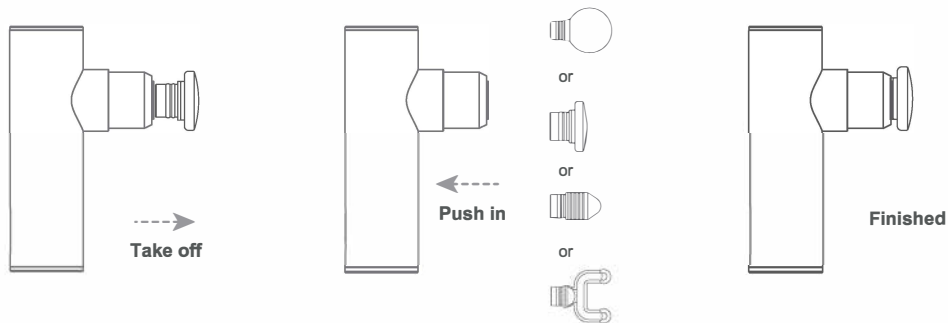Turn on/off


- Press and hold the button for 3 seconds to turn on the device.
- Press and hold the button for 3 seconds to turn off the device.

Speed Adjustment

- Level 1 will be on by default when the device is turned on.
- Each press will increase the speed to a higher level.
- When the speed reaches Level 3, press once more and the speed will return to Level 1.

Attention:

The device is equipped with 10 minutes of working timer to avoid damage to the muscle and fascia caused by excessive massage.

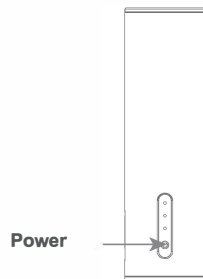Indicator and Blinker

- The speed level indicator shows the vibration speed.
- Level 1 will be on by default when the device is turned on.
- The speed will be changed by each short-press sequentially and cyclically.


3 Levels of Speed Intensity

Level 1	Level 2	Level 3
1800r/min	2600r/min	3200r/min

Speed Level Indicator

Blinker/Power Indicator

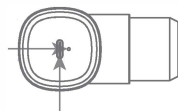
Charging Port

The red light is on when the device is charging.

The white light is on when the device finishes charging.

About Charging

- Fully charging the device before the first time you use it.
- The red light will be on when the device is charging, and the white light will be on when the charging is finished.
- Full charging will take around 3h, and the run time will last around 8h when the level is at the highest.
- When the battery level is less than 10%, the red light will flash and please charge it as soon as possible.

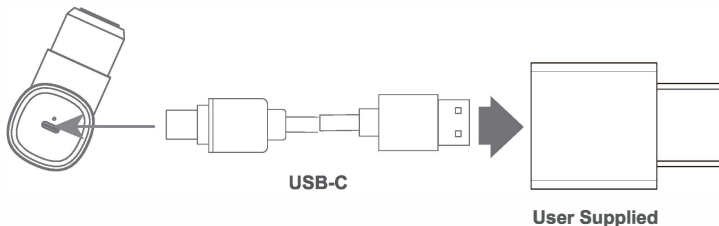Note: The charger should be 5V/2A.

*For your security, do not use the device when it is charging.

*For a long service life, please charge the device as soon as possible once the red light flashes.

Contraindication

PLEASE DO NOT USE THE DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY.

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health.

Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above.

Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

Warning & Caution

PLEASE DO NOT USE THE DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY.

- DO NOT immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- DO NOT remove screws or attempt to disassemble.
- DO NOT operate continuously for more than one hour. Allow device to rest for 30 minutes before reuse.
- Unplug the charger from device battery after charging or prior to use.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

Safety Instructions

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE AND PERSONAL INJURY, THIS PRODUCT MUSTBE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.

- For adult use only.
- Please keep this device away from children or store it in places where children cannot reach.
- This device is not food. Please do not put the head attachment or any parts into your mouth.

- Use only on the body's soft tissue as desired. Do not use on head or any hard or bony area of the body.
- Do not apply this device to the following body parts: sex organs, face, eyes, injuries after surgery, any kind of infection part, infected or damaged skin, injured bones and joints.
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- In case of pain or other discomfort, please stop using immediately.
- Use only the head attachment which Hi5 provided and use on applicable body area.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Please keep fingers, hair or any other body part away from the shaft and backside of the head attachment as pinching may occur.
- Do not place any objects into the ventilation holes of the device.
- Carefully examine the device and battery before each use.
- Do not tamper or alter the device.
- Never leave the Vibration Massage Device working or charging unattended.

Battery Disposal

Before discarding the device, please use a suitable tool to open the metal part of the bottom, and then pick up the electronic board assembly so that you can take out the battery inside to discard (the user cannot replace the battery by himself):

- The battery must be removed from the device before discarding the device.
- When removing the battery, make sure the device is powered off.
- The battery should be disposed of according to the national garbage classification and hazardous garbage classification.

Warranty

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OF WORKMANSHIP, HI5 WILL REPAIR OR REPLACE THE PRODUCT, OR SUPPLY NECESSARY COMPONENTS FREE OF CHARGE.

This Warranty Excludes:

- a) Damage caused by accident, abuse, mishandling, or transport;
- b) Units subjected to unauthorized repair;
- c) Improper use of this device;
- d) Damage exceeding the cost of the product;

- e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises;
- f) Failure to provide the dated proof of purchase.

● Packing List